Resin-based Dental Restorative Material

# OMNICHROMA BLOCKER

Flow

CE 1639

SYRINGE

## Instructions for Use

### Tokuyama Dental Corporation

38-9, Taitou 1-chome, Taitou-ku, Tokyo 110-0016, Japan

Tel: +81-3-3835-7201

URL: <https://tokuyama-dental.com>

|          |    |             |    |             |    |            |    |
|----------|----|-------------|----|-------------|----|------------|----|
| ENGLISH  | 2  | БЪЛГАРСКИ   | 5  | ČESKY       | 10 | DEUTSCH    | 14 |
| DANSK    | 18 | EESTI KEEL  | 22 | ESPAÑOL     | 26 | SUOMI      | 30 |
| FRANÇAIS | 34 | ΕΛΛΗΝΙΚΑ    | 38 | HRVATSKI    | 43 | MAGYAR     | 46 |
| ITALIANO | 51 | LIETUVIŲ K. | 55 | LATVISKI    | 59 | NEDERLANDS | 63 |
| NORSK    | 67 | POLSKI      | 71 | PORTUGUÊS   | 75 | ROMÂNĂ     | 79 |
| SVENSKA  | 83 | SLOVENČINA  | 87 | SLOVENŠČINA | 91 | TÜRKÇE     | 95 |

This product is certified as a medical device under the legislation of the European Union on medical devices by SGS CE1639, exclusively for the indication(s) shown in this IFU. Other non-medical uses ascribed to this product are not within the scope of CE certification, and users should be aware product performance and/or safety has not been evaluated by SGS for those purposes.

Carefully read all information, precautions, notes and materials before using.

## **PRODUCT DESCRIPTION AND GENERAL INFORMATION**

- 1) OMNICHROMA BLOCKER FLOW is a low viscosity, light-cured, radiopaque composite resin for use in anterior and posterior restorations. Direct placement of OMNICHROMA BLOCKER FLOW into a prepared cavity allows for easy handling. It has radiopacity equal to or greater than 1mm aluminum, and is classified as a Type 1 and Class 2 (Group 1) material by ISO 4049.
- 2) OMNICHROMA BLOCKER FLOW contains 71% by weight (57% by volume) of spherical silica-zirconia filler and composite filler. The spherical filler contained in OMNICHROMA BLOCKER FLOW (mean particle size: 0.2  $\mu\text{m}$ , particle size range: 0.1 to 0.3  $\mu\text{m}$ ) facilitates excellent gloss retention and wear resistance.
- 3) OMNICHROMA BLOCKER FLOW contains Bisphenol A di(2-hydroxy propoxy) dimethacrylate (Bis- GMA), Bisphenol A polyethoxy methacrylate (Bis-MPEPP), Triethylene glycol dimethacrylate (TEGDMA), 1,6-bis(methacryl-ethyloxycarbonylamino)trimethyl hexane (UDMA), Mequinol, Dibutyl hydroxyl toluene and UV absorber.
- 4) OMNICHROMA BLOCKER FLOW is available in SYRINGE.
- 5) OMNICHROMA BLOCKER FLOW is designed for masking slight stain and discoloration. OMNICHROMA BLOCKER FLOW can be used as a thin layer at the lingual cavity wall of extensive class III and IV restorations or reconstruct a highly opaque tooth.

## **INDICATIONS**

For use as a tooth shade resin in dental procedures, such as:

- Direct anterior and posterior restorations
- Cavity base or liner
- Repair of porcelain/composite

## **CONTRAINDICATIONS**

OMNICHROMA BLOCKER FLOW contains methacrylic monomers and UV absorber. DO NOT use OMNICHROMA BLOCKER FLOW for patients allergic to or hypersensitive to methacrylic and related monomers, UV absorber or any of the other ingredients.

## **PRECAUTIONS**

- 1) DO NOT use OMNICHROMA BLOCKER FLOW for any purpose other than those listed in these instructions. Use OMNICHROMA BLOCKER FLOW only as directed herein.
- 2) OMNICHROMA BLOCKER FLOW is designed for sale and use by licensed dental care professionals only. It is not designed for sale nor use by non-dental care professionals.
- 3) DO NOT use OMNICHROMA BLOCKER FLOW if the safety seals are broken or appear to have been tampered with.
- 4) If OMNICHROMA BLOCKER FLOW causes an allergic reaction or hypersensitivity, discontinue its use immediately.
- 5) Use examination gloves (plastic, vinyl or latex) at all times when handling OMNICHROMA BLOCKER FLOW to avoid the possibility of allergic reactions from methacrylic monomers or UV absorber.

Note: Certain substance/materials may penetrate through examination gloves. If OMNICHROMA BLOCKER FLOW comes in contact with the examination gloves, remove and dispose of the gloves, and wash hands thoroughly with water as soon as possible.

- 6) Avoid contact of OMNICHROMA BLOCKER FLOW with eyes, mucosal membrane, skin and clothing.
  - If OMNICHROMA BLOCKER FLOW comes in contact with the eyes, thoroughly flush

eyes with water and immediately contact an ophthalmologist.

- If OMNICHROMA BLOCKER FLOW comes in contact with the mucosal membrane, wipe the affected area immediately, and thoroughly flush with water after the restoration is completed.
  - If OMNICHROMA BLOCKER FLOW comes into contact with the skin or clothing, immediately saturate the area with an alcohol soaked cotton swab or gauze.
  - Instruct the patient to rinse his mouth immediately after treatment.
- 7) OMNICHROMA BLOCKER FLOW should not be ingested or aspirated. Ingestion or aspiration may cause serious injury.
  - 8) To avoid the unintentional ingestion of OMNICHROMA BLOCKER FLOW, do not leave it unsupervised within the reach of patients and children.
  - 9) Clean the placement instruments and brushes with alcohol after use.
  - 10) When using a light-curing unit, protective eye shields, glasses or goggles should be worn at all times.
  - 11) This product is designed to be used at room temperature (18 - 30°C / 62 - 84°F). Allow product to reach room temperature prior to use. Cold material may be difficult to extrude. After removing from the refrigerator, allow the product to sit 15 minutes or until it reaches room temperature.
  - 12) Be aware that when a restoration is chipped from misaligned occlusion or bruxism (clenching, grinding or tapping), the repaired restoration may chip again.

## ■ INDICATIONS FOR FILLING AND CURING

OMNICHROMA BLOCKER FLOW is designed to be cured by either a halogen or LED curing-light with a wavelength of 400-500 nm. Be sure to light-cure OMNICHROMA BLOCKER FLOW extra-orally and check the time needed for complete hardening of OMNICHROMA BLOCKER FLOW with your light-curing unit before performing the bonding procedure. The following table summarizes the relationship between curing time and increment depth.

Relationship between curing time and increment depth:

| Light type | Intensity<br>(mW/cm <sup>2</sup> ) | Curing time<br>(seconds) | Increment depth<br>(mm) <sup>1)</sup> |
|------------|------------------------------------|--------------------------|---------------------------------------|
| Halogen    | 400                                | 20                       | 1.5                                   |
|            |                                    | 40                       | 2.0                                   |
|            | 600                                | 20                       | 1.6                                   |
|            |                                    | 10                       | 1.5                                   |
|            |                                    | 20                       | 1.7                                   |
| LED        | 400                                | 20                       | 1.5                                   |
|            |                                    | 40                       | 1.9                                   |
|            | 600                                | 20                       | 1.4                                   |
|            |                                    | 10                       | 1.3                                   |
|            | 800                                | 20                       | 1.6                                   |
|            |                                    | 10                       | 1.5                                   |
|            | 2300                               | 6                        | 1.9                                   |

- 1) Increment depth was determined on the basis of test results performed in accordance with “depth of cure” of ISO 4049.

## ■ STORAGE

- 1) Store OMNICHROMA BLOCKER FLOW at temperatures between 0 - 25°C (32 - 77°F).
- 2) AVOID direct exposure to light and heat.
- 3) DO NOT use OMNICHROMA BLOCKER FLOW after the indicated date of expiration on

the syringe package.

## ■ DISPOSAL

To safely dispose of excess OMNICHROMA BLOCKER FLOW, extrude unused portion from SYRINGE and light-cure before disposal.

Do not allow to enter drains, sewers or watercourses.

Must not be disposed of together with household garbage.

Follow local instructions for disposal.

Refer to the packaging waste disposal guide shown on the end of this IFU.

## ■ CLINICAL PROCEDURE

### 1. Cleaning

Thoroughly clean the tooth surface with a rubber cup and a fluoride-free paste then rinse with water.

### 2. Shade Selection

Select OMNICHROMA BLOCKER FLOW as the initial layer under OMNICHROMA or OMNICHROMA FLOW, when masking slight stain or extensive Class III or IV direct cavity restorations.

### 3. Isolation

A rubber dam is the preferred method of isolation.

### 4. Cavity Preparation

Prepare the cavity. In case where no cavity preparation has been made (caries-free cervical defects), clean the tooth surface with a rubber cup and a fluoride-free cleaning paste. Rinse thoroughly with water.

- In the case of porcelain/composite repairs, roughen the surface with a bur or a diamond point to prepare the area for adhesion; apply phosphoric acid etch for cleaning; rinse thoroughly with water; air dry thoroughly and treat with a silane coupling agent, followed by a metal primer if metal is exposed on the fractured surfaces according to manufacturer's instructions.

### 5. Pulp Protection

Calcium hydroxide should be applied if the cavity is in close proximity to the pulp. DO NOT USE EUGENOL-BASED MATERIALS to protect the pulp as these materials will inhibit curing OMNICHROMA BLOCKER FLOW. OMNICHROMA BLOCKER FLOW is not indicated for direct pulp capping.


### 6. Bonding System

Apply a bonding system according to its manufacturer's instructions.

- If using self-cured or dual-cured bonding systems, please confirm that such bonding systems are compatible with the light-cured composite resins being used.
- In case of porcelain/composite repairs, be sure to condition the surface of porcelain/composite with a silane coupling agent BEFORE performing the bonding procedure.

### 7. Dispensing

- Hold the nozzle and turn the cap counterclockwise for removal. [1]
- Hold the nozzle and place the Dispensing Tip (metal tip) onto it. Turn the Dispensing Tip clockwise until it is securely locked. [2]

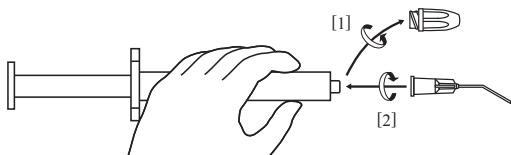- Accessories and Dispensing Tips are exclusively designed for the SYRINGE of OMNICHROMA BLOCKER FLOW. Use the attached Dispensing Tips (metal tips).
- Dispensing Tips (TOKUYAMA Dispensing Tip) are also sold separately.
- Use uncontaminated gloves when handling the syringe and attaching the Dispensing Tip.

## **8. Filling and contouring**

Fill the cavity with OMNICHROMA BLOCKER FLOW as an initial layer. Cavities can be directly filled from the Dispensing Tip, or indirectly filled with an instrument after extruding the paste on the mixing pad. Increments should not exceed the indicated curing depth at a time (please refer to the aforementioned table).

- DO NOT apply excessive force to extrude the material from the syringe immediately after removal from the refrigerator.
- DO NOT mix OMNICHROMA BLOCKER FLOW with other brands of resin composite to avoid incomplete cure, entrapment of air bubbles or losing its color adjustment.
- To avoid cross infection, after extruding the paste, remove and discard the Dispensing Tip, wipe the nozzle with gauze and replace the cap immediately.
- This product is provided non-sterile and is not intended to be sterilized by the user.

## **9. Curing**

Light-cure each increment for at least the indicated time (see aforementioned table).

- Follow the filling and curing instructions of OMNICHROMA FLOW/ OMNICHROMA to be layered over the cured OMNICHROMA BLOCKER FLOW.

## **10. Finishing**

Shape and polish the restoration according to instructions for OMNICHROMA/ OMNICHROMA FLOW used as the surface layer of the restoration.

The user and/or patient should report any serious incident that has occurred in relation to the device to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.

**IMPORTANT NOTE:** The manufacturer is not responsible for damage or injury caused by improper use of this product. It is the personal responsibility of the dental professional to ensure the product is suitable for application before use.

Specifications are subject to change without notice. When the product specification changes, the instructions and precautions may change also.

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# Packaging waste disposal guide

(According to European Commission Decision 97/129/EC)

## ENGLISH

### Packaging waste disposal guide

|                           |                                                                                          |
|---------------------------|------------------------------------------------------------------------------------------|
| SYRINGE (cap and body)    | Dispose of as industrial waste.                                                          |
| Packaging bag for tip     | LDPE4 (Low density polyethylene):<br>Recycling. Follow local instructions for disposal.  |
| Paperbox and/or partation | PAP21 (Non-corrugated fibreboard):<br>Recycling. Follow local instructions for disposal. |
| IFU and/or guide          | PAP22 (Paper):<br>Recycling. Follow local instructions for disposal.                     |

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