


**AUTOMATIC  
EFFICIENCY**

**SPEED, PRECISION,  
RELIABILITY IN EVERY  
APPLICATION.**


**ESTELITE  
 $\Sigma$ QUICK**

**Maximum modelling freedom,  
minimum polymerisation time.**

With Estelite Sigma Quick, achieving aesthetic and durable restorations is quick and easy. Thanks to its innovative technology, this composite offers superior handling and a long-lasting shine, reducing times without compromising on quality. **Efficiency you can count on.**

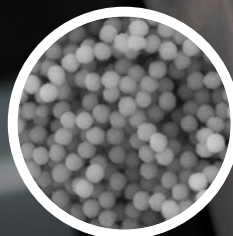

ESTELITE SIGMA  
20,000X

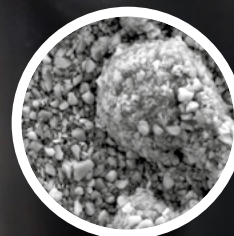

CERAM.X SPHERETEC ONE  
UNIVERSAL\* 5,000X

\* It is not a registered trademark  
Tokuyama Dental Corporation

# ESTELITE Σ QUICK


Authentic  
smiles,  
always





## Features

- Ultra-fast curing time  
3 seconds with plasma lamp  
(maximum 2 mm thickness)
- Extended working time  
90 seconds in 10,000 lx light  
environment
- Easy to polish
- Long-lasting brightness
- True chameleon effect
- High abrasion resistance
- Extremely low polymerization  
shrinkage (1.3%)
- Radiopaque
- Anterior and posterior  
restorations
- Closure of diastemas
- Masking of discolorations
- Veneer in composite
- Cementation of indirect  
restorations
- Characterizations and color  
customizations for direct and  
indirect resin-based restorations

Dr.  
Bilal Yaşa

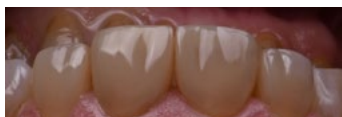

Dr.  
Nicolò Barbera

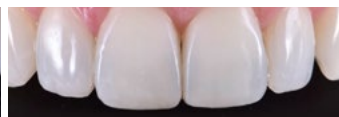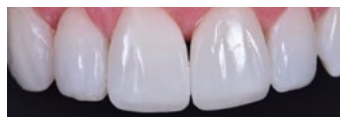

Winning technology for  
more than 10 years.

**DA**  
DENTAL  
ADVISOR



2010 / 2024 TOP AWARD COMPOSITE



## Available packaging

# Estelite Sigma Quick Intro Kit

- 6 Estelite Sigma Quick 3.8g syringes (colors: A1, A2, A3, A3.5, OA2, OPA2)
- 1 Bond Force II Bottle 1ml 50 applicators
- 1 mixing cup Shade Guide



# Estelite Sigma Quick 3 Syringe Kit

- 3 Estelite Sigma Quick syringes 3.8g (colors: A2, A3, OA3)

## Spare parts

- 1 Estelite Sigma Quick 3.8g syringe
- 20 capsules (PLT) Estelite Sigma Quick, 0.20g Estelite Sigma Quick Shade Guide

## Color classification

- Enamels



- Body



- Effect masses



- Dentines



## Tips and tricks

Use OPA2 in classes III and IV as the first layer.


Did you know that OPA2 can be used also with Estelite Asteria?

# Sol-Gel Method


SINCE 1981

The Supra-Nano sized spherical filler for simple, long-lasting and wonderful results.

COPYING IS EASY, INNOVATING IS TOKUYAMA.

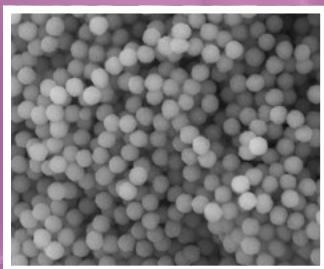

ESTELITE  
SIGMA QUICK  
20,000X

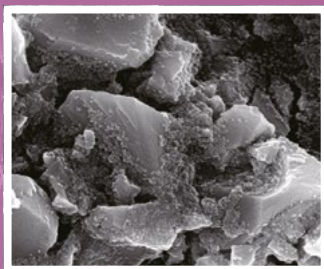

GÆNIAL «GC»\*  
20,000X

\* Not a registered  
trademark of  
Tokuyama Dental  
Corporation

Tokuyama, thanks to the patented Sol-Gel Method, is able to determine in the laboratory the shape and size of its  $\text{SiO}_2$  and  $\text{ZrO}$ -based fillers <sup>2</sup> (silica and zirconium oxide), resulting in a truly spherical filler called Supra Nano Spherical Filler.

It is a compound made thanks to a chemical reaction of metal alkoxides. Thanks to this method, the Supra-Nano molecular dimensions give the material characteristics that allow restorations with very smooth surfaces and maximum final brilliance.



AN EXCLUSIVE PATENT  
TOKUYAMA DENTAL SINCE 1981

 **Tokuyama**  
Innovating Tomorrow's Dentistry, Today.

Tokuyama Dental Italy srl  
tel. +39 0444 659650  
fax +39 0444 750345  
info@tokuyama.it



[www.tokuyama.it](http://www.tokuyama.it)

DPSQEN01