

AUTOMATIC EFFICIENCY

SPEED, PRECISION,
RELIABILITY IN EVERY
APPLICATION.


ESTELITE Σ QUICK

**Maximum modelling freedom,
minimum polymerisation time.**

With Estelite Sigma Quick, achieving aesthetic and durable restorations is quick and easy. Thanks to its innovative technology, this composite offers superior handling and a long-lasting shine, reducing times without compromising on quality. **Efficiency you can count on.**



ESTELITE ASTERIA
20,000X

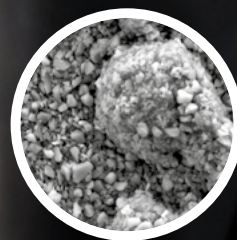

CERAM.X SPHERTEC ONE
UNIVERSAL* 5,000X

ESTELITE Σ QUICK


Authentic
smiles,
always





Features

- Ultra-fast curing time
3 seconds with plasma lamp
(maximum 2 mm thickness)
- Extended working time
90 seconds in 10,000 lx light
environment
- Easy to polish
- Long-lasting brightness
- True chameleon effect
- High abrasion resistance
- Extremely low polymerization
shrinkage (1.3%)
- Radiopaque
- Anterior and posterior
restorations
- Closure of diastemas
- Masking of discolorations
- Veneer in composite
- Cementation of indirect
restorations
- Characterizations and color
customizations for direct and
indirect resin-based restorations

Dr.
Bilal Yaşa

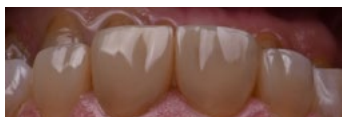

Dr.
Nicolò Barbera

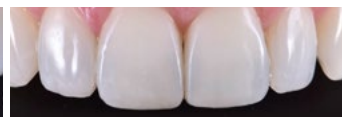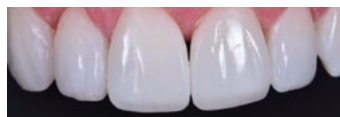

Winning technology for
more than 10 years.

DA
DENTAL
ADVISOR



2010 / 2024 TOP AWARD COMPOSITE



Available packaging

Estelite Sigma Quick Intro Kit

- 6 Estelite Sigma Quick 3.8g syringes (colors: A1, A2, A3, A3.5, OA2, OPA2)
- 1 Bond Force II Bottle 1ml 50 applicators
- 1 mixing cup Shade Guide



Estelite Sigma Quick 3 Syringe Kit

- 3 Estelite Sigma Quick syringes 3.8g (colors: A2, A3, OA3)

Spare parts

- 1 Estelite Sigma Quick 3.8g syringe
- 20 capsules (PLT) Estelite Sigma Quick, 0.20g Estelite Sigma Quick Shade Guide

Color classification

- Enamels



- Body



- Effect masses



- Dentines



Tips and tricks

Use OPA2 in classes III and IV as the first layer.


Did you know that OPA2 can be used also with Estelite Asteria?

Sol-Gel Method


SINCE 1981

The Supra-Nano sized spherical filler for simple, long-lasting and wonderful results.

COPYING IS EASY, INNOVATING IS TOKUYAMA.

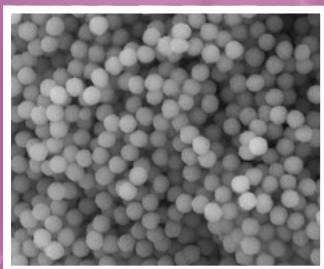

ESTELITE
SIGMA QUICK
20,000X

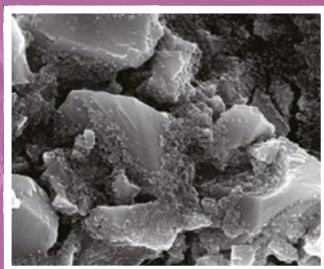

GÆNIAL «GC»*
20,000X

* Not a registered
trademark of
Tokuyama Dental
Corporation

Tokuyama, thanks to the patented Sol-Gel Method, is able to determine in the laboratory the shape and size of its SiO_2 and ZrO -based fillers ² (silica and zirconium oxide), resulting in a truly spherical filler called Supra Nano Spherical Filler.

It is a compound made thanks to a chemical reaction of metal alkoxides. Thanks to this method, the Supra-Nano molecular dimensions give the material characteristics that allow restorations with very smooth surfaces and maximum final brilliance.



AN EXCLUSIVE PATENT
TOKUYAMA DENTAL SINCE 1981

 **Tokuyama**
Innovating Tomorrow's Dentistry, Today.

Tokuyama Dental Italy srl
tel. +39 0444 659650
fax +39 0444 750345
info@tokuyama.it



www.tokuyama.it

DPSQEN01